EN: This is the simplified manual. You can consult the detailed manual via our official website or by clicking on the link at the top of all our emails.

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ACUPEO®

MANUAL

EN: Remember, positive feedback is essential to help us continue to make a living from what we love and to improve our products and services in the future. Thank you!

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ACUPEO® PULSE TENS MASSAGER

Rechargeable Version

1. Remove the plastic protection from the pad.
2. Place the device on the area of the body to be massaged, avoiding the heart, brain, eyes, and genital areas.
3. Turn on the device by pressing and holding the "ON-INC" button.
4. Select a mode by pressing the "PROGRAM" button. You can choose from 6 modes.
Acupuncture: OXXX Kneading: X0XX Tapping: 00XX
Physiotherapy: XX0X Beat: XXX0 Stimulation: XX00.
5. Choose from 10 intensity levels using the "ON-INC" button to increase and "OFF-DEC" to decrease.
6. To turn off the device, simply press and hold the "OFF-DEC" button.
7. At the end of the session, clean the pad with water and put back its plastic protection to extend its lifespan.

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
→If you have any problems getting the device to work, you can contact us via our official website or by clicking on the contact link at the top of all our emails.

→If you would like to order other pads compatible with ACUPEO®, visit the menu on our official website or click on the order link at the top of all our emails.

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ACUPEO® all right reserved

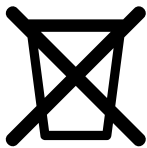Our TENS emits a small current: Do not use the device on the heart, eyes, genitals, brain or throat; Do not use the device in the shower or bath; Recent pacemakers are often compatible but you must ask your cardiologist, as is the case for other heart devices/problems; You must ask the doctor's advice if you are pregnant or if you have a particular medical condition such as epilepsy.

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